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## TM-5704

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**MASSAGING FOOT WARMER**

**MODEL: TM-5704**

inside the box and either store or dispose of safely. Make sure you save these instructions for future reference.

## **IN THE BOX**

- Foot warmer Massage unit
- Controller
- Power adaptor
- Instruction manual

## **FEATURES**

- Hand held remote control
- Heat setting
- Variable massage control
- Comfortable soft fleece lining

## **SAFETY INSTRUCTIONS**

- When using electrical products, basic safety precautions should always be followed.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities or lack of experience or knowledge. They should be given supervision and instruction in the use of the appliance by a person responsible for their safety.
- Check that the voltage indicated on the rating plate corresponds with that of the local network before connecting the product to the mains supply.
- Do not attempt to make any modifications in any way to the product.
- Ensure the product is not too close to any heat emitting appliances.
- Do not leave the appliance unattended for a long period of time when switched on.
- Do not operate this product under a blanket or a pillow.
- Do not immerse the product in water or any other liquids.
- Repairs to electrical appliances should only be performed by a qualified electrician.

- Improper repairs may place the user at serious risk.
- This appliance is intended for domestic use only.
- Do not operate this product after a malfunction, or after being dropped or damaged in any way. If the adaptor becomes damaged it must be replaced by a qualified service engineer.
- Do not operate where aerosol products are being used or where oxygen is being administered.
- Do not use this appliance for anything other than its intended use.
- Do not use the appliance outdoors.
- Switch off this product when not in use.
- This item is NOT machine washable.


Note to User:

- Do not use the massager on the following body areas in these instances:
  - Any area of the body that is swollen, burnt or inflamed or where skin eruptions or ulcerated sores are present.
  - Any area of the body that is anaesthetised unless approved by your doctor.
  - Legs with varicose veins.
  - Calves or legs with unexplained pain.
  - Do not use for more than 15 minutes at a time.

Massage should be pleasant and comfortable. Should pain or discomfort occur, discontinue use immediately and/or consult your doctor.

## DESCRIPTION OF PARTS
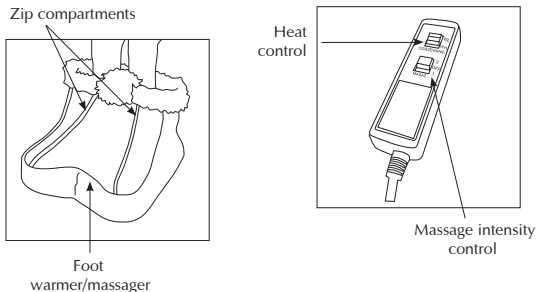
## INSTRUCTIONS FOR USE

Now you can be warm and cosy even on the coldest of nights, with the Massaging Foot Warmer. Settle sore or chilly feet into the soft, fleece lining and let your cares melt away. The remote control allows adjustment of the heat and intensity of the massage as desired; you can even choose to use the heat or massage functions separately.

## HOW TO USE

- Place the foot massager on a flat surface.
- Plug the power adaptor into a suitable socket outlet and connect the receptacle end into the power input of the foot massager. Switch on at the socket outlet if switchable.
- Place both feet inside the boot for warmth.
- For additional warmth, turn the heat setting on the remote to the On position.
- Note: Do not use the heat feature excessively in any one session.
- To activate the vibrating massage, set the massage setting to On.
- Select the intensity of the massage from Low to High using the remote control.

Note: Do not use the massage feature for more than 15 minutes in any one session.

## TECHNICAL SPECIFICATION

- Electrical Supply
- Foot Warmer: Input 12V DC
- Adaptor: Input 230V SOHZ Output 12V DC 1 000mA

## TROUBLESHOOTING

Problem	Cause	Solution
No heat from appliance.	Adaptor not fitted cor-rectly.	Ensure power supply is securely connected.
	Foot warmer unzipped.	Ensure all zips are adequately fastened.
Vibrating function not active.	Adaptor not correctly fitted.	Ensure power supply is securely connected.

## CARE & MAINTENANCE

- Always unplug the appliance from the mains power supply before cleaning the unit.
- Never immerse the product in water or any other liquids.
- Clean the outer surface of the appliance with a soft,damp cloth and a mild washing liquid if necessary. Do not use abrasive cleaning fluids as these could damage the outer surface.
- Never machine wash this product.
- Allow all parts to dry thoroughly before connecting to the power supply.

The crossed out wheelie bin symbol on this item indicates that this appliance needs to be disposed of in an environmentally friendly way when it becomes of no further use or has worn out. Contact your local authority for details of where to take the item for re-cycling.

## Disposal

Batteries and electronic instruments must be disposed of in accordance with the locally applicable regulations, not with domestic waste.

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